

Daisy

Versatile, ergonomic, and secure storage



 Lithium-ion battery

 Electric height adjustment

This cart is designed to be used with different computers and devices. Our medical-grade Infinity all-in-one PC is shown on the cart above and features:

- A large, intuitive touchscreen.
- Dual hot-swappable batteries for 24/7 constant power.
- Powerful CPU, ultra-fast SSD, and up to 64GB system memory.
- Unique battery management software with remote monitoring.
- Front panel is IP65 rated for protection from water and dust.
- 3-year warranty.

For more info on our Infinity all-in-one PCs or to enquire about a cart bundle package, please get in touch with our sales team, who will be happy to help. Our contact details are below.

Features

Combining practical function with an infection control-friendly design, Daisy is ideal for light medication rounds or for keeping essential supplies to hand during clinical rounds.

- ✓ Electric height adjustment to change height at the touch of a button, with independent height adjustment for your monitor.
- ✓ A slide-out side tray offers additional medication preparation space, while a lockable metal drawer keeps supplies organised and secure. Digi, keyed, and RFID lock options are available.
- ✓ Modular design allows for the addition of a range of accessories.
- ✓ Retractable keyboard conveniently tucks out of the way during medication prep and dispensing.
- ✓ Work surface is durable and resistant to hospital cleaning agents.

Benefits

- Small footprint makes the cart easy to position in busy bays.
- Large work surface and slide out tray provide ample medication preparation space.
- Conveniently placed grab handles at the side of the cart provide an unobscured view when moving the cart to a new area.
- The wider work surface design offers protection from damage for your screen while being moved to a new area.
- LED battery gauge at front of cart provides an at-a-glance view of battery power status from across the room.
- Easily accessible PC enclosure and cable management.
- Up to two optional side rails provide versatility for mounting a range of accessories, such as label printers, barcode scanners, and sharps bins.
- Designed for infection control and ease of cleaning.
- Comprehensive 5-year cart warranty.

Areas of use

Daisy is ideal for providing access to electronic patient records as well as secure access to medication and supplies in areas such as:

- A&E / ED
- ICU
- Theatres
- Wards

Compatible with



Infinity
all-in-one PC



Micro
PC




Laptop


Daisy

Versatile, ergonomic, and secure storage

Drawings

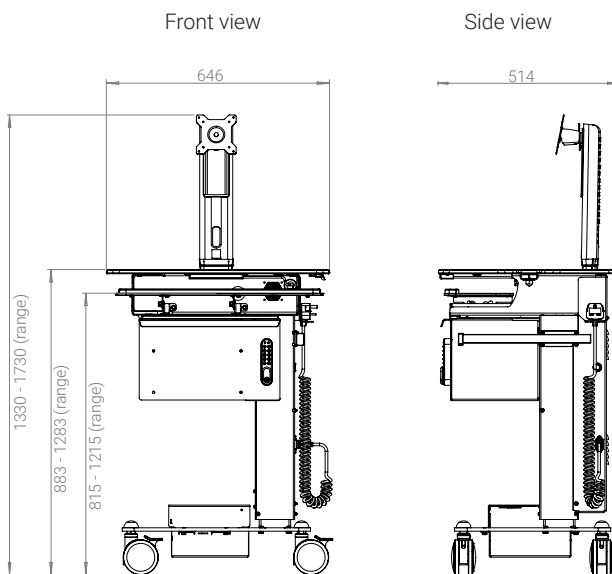

PC caddy tray location and dimensions

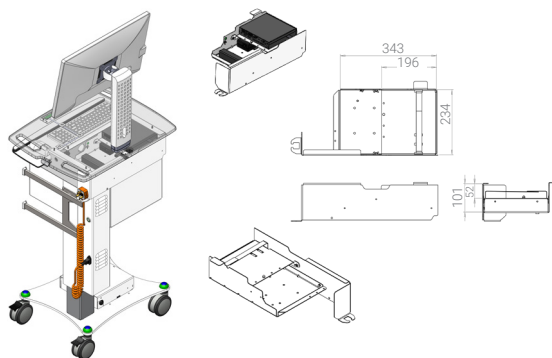

Main specifications

Work surface height (extended)	1283 mm
Work surface height (compressed)	883 mm
Work surface dimensions	646 mm (width) x 514 mm (depth)
Storage drawer internal dimensions	277 mm (width) x 380 (length) x 138 mm (depth)
Weight (cart base only)	63 kg
Weight (max weight possible)	78 kg
PC caddy dimensions	386 mm (width) x 234 mm (depth) x 103 mm (height)
PC shelf dimensions	198 mm (width) x 234 (depth) x 49 mm (height)

Approved cleaning solutions

- Quaternary ammonium chloride compounds (such as AirX109)
- Bleach-Household strength (1:10 minimum dilution strength)
- Ethyl or Isopropyl alcohol-70%
- Hydrogen Peroxide
- Phenolic disinfectants (such as Lysol, WexCide)
- Ovation
- Fantastik
- Formula 409
- 1.6% Aqueous Ammonia
- "Green Soap" USP (United States Pharmacopoeia)
- Clinell Wipes
- Chlor Clean



Direct purchase



Made in the UK



Project management



Device configuration



Phased delivery



Instant support

Common accessories



Printer



Scanner



Medical keyboard & mouse

Version 1.0

